


Catheter Accessories

Titanium Catheter Adapter 2-Part Order No. TA-200




Catheter Plastic Adapter
Order No. PCA-1010

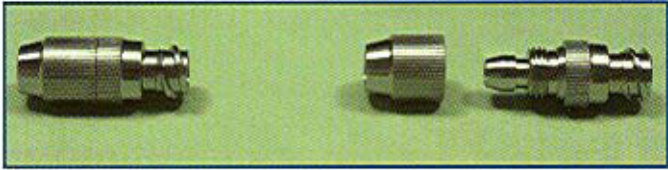


Titanium Catheter Extender
Order No. TE-1800 (Adult) TE-1800P (Pediatric)



Catheter Clamp
Order No. PCA-1030

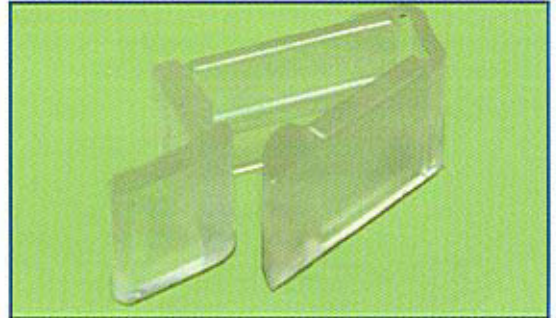


Trocar (Faller Stylet) Order No. FS-401
Autoclavable (designed to fit Tenckhoff catheters)

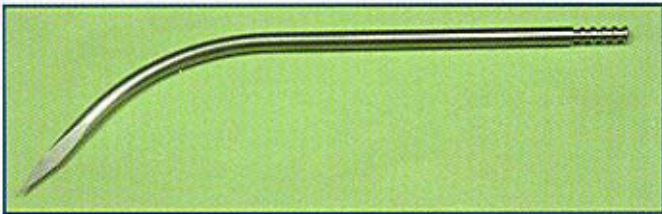

Injection Cap
Order No. ICL-100



Stylet (Autoclavable)
Order No. SST - XXX (length in mm)

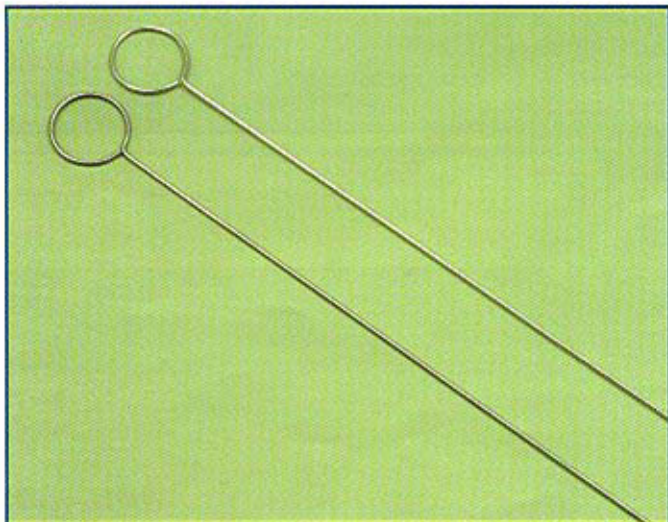


QC Cap (Prefilled with Povidone Iodine)
Order No. QC-200
(used to cap off PD catheters in between CAPD exchanges) (helps to reduce infection)




Catheter Accessories con't.


Trocar (Autoclavable) Order No. TR-400
(designed to fit standard Tenckhoff
or Reinforced Tenckhoff catheters)

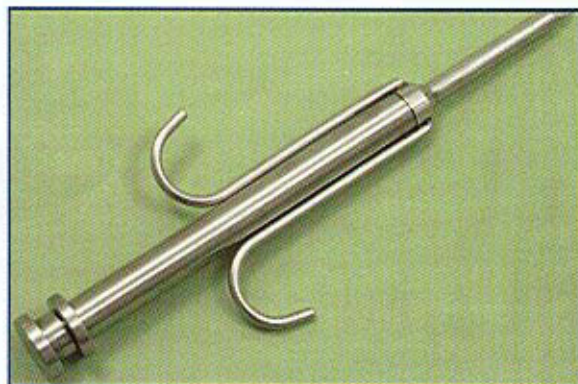

Pull-Apart Introducer Set Order No. PS-16FR
(used to insert catheter percutaneously)

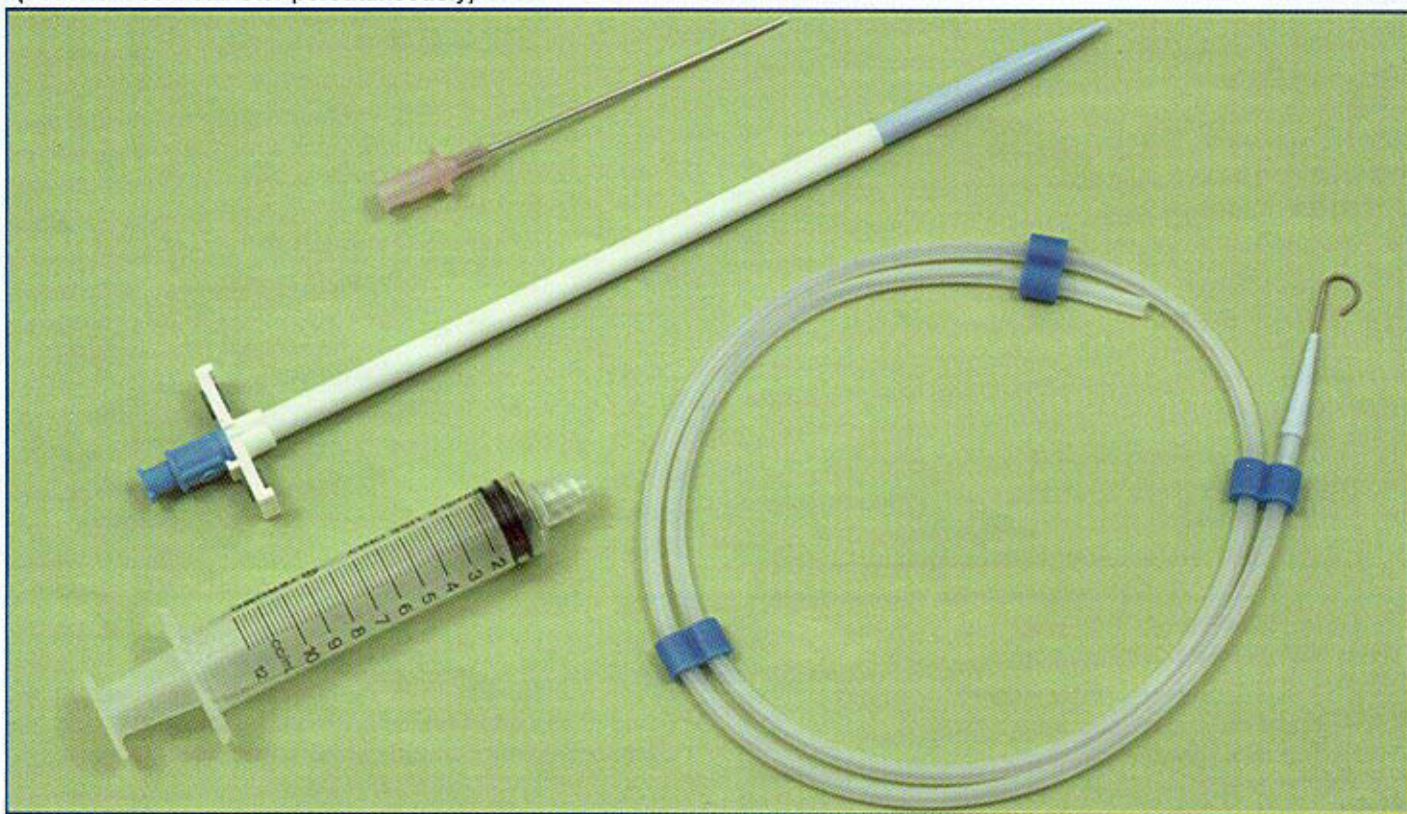

Peritoneal Tunneling Stylet (Autoclavable)
Order No. PS-423
(designed to fit Tenckhoff Catheters)

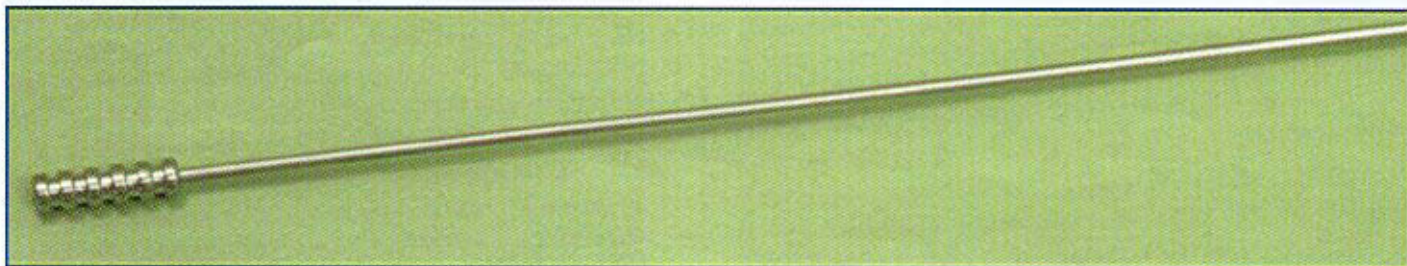

Catheter Accessories con't.

GUIDEWIRE

Order No. GW-3870

J/Straight

0.97mm x 70cm


(0.038" x 27.6")

Order No. GW-3570

J/Straight

0.89mm x 70cm

(0.035" x 27.6")

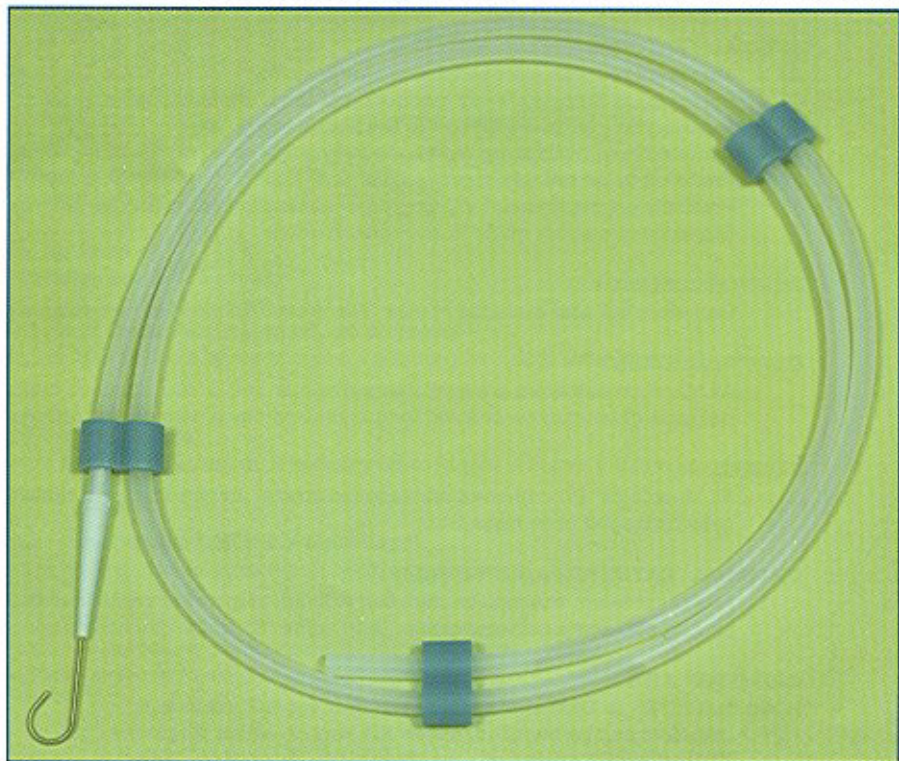

INTRODUCER NEEDLE

Order No. CN-18G7 (18 Ga x 7cm)


Order No. CN-16G7 (16 Ga x 7cm)



DILATOR

Order No. DL-10F6 (10 Fr x 6")

Order No. DL-12F8 (12 Fr x 8")

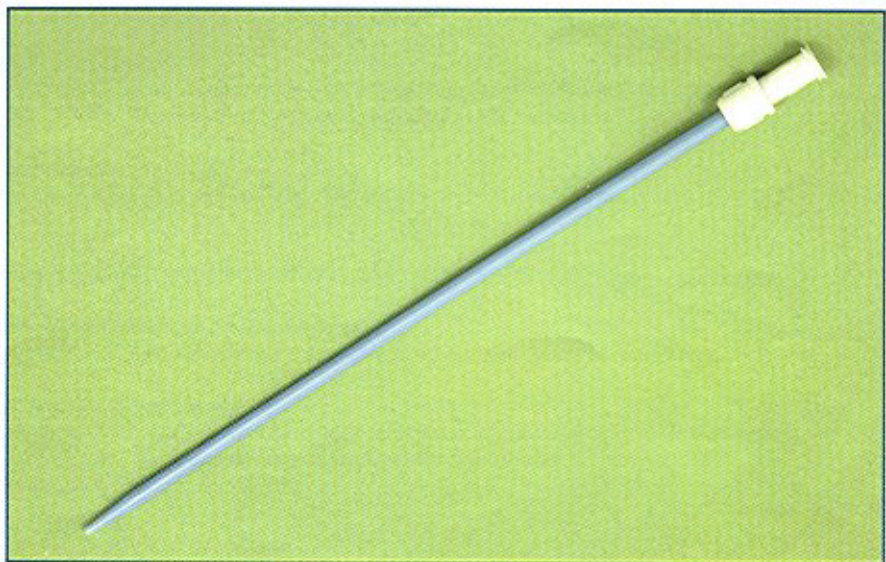

Peritoneal Dialysis Catheter Accessories

Catheter Adapter PCA-1010 (Sterile, individually packed, 30 pkgs/cs)

... is an inexpensive alternative to the Titanium Adapter. It allows peritoneal catheters to connect securely to CAPD administration sets.

Injection Cap ICL-100 (Sterile, individually packed, 30 pkgs/cs)

... used to cap off the catheter adapter. It is self sealing and provides for administering medication.

Catheter Clamp PCA-1030 (Sterile, individually packed, 30 pkgs/cs)

... allows easy clamping of the peritoneal catheter.

Catheter Adapter Kit PCA-1000K (Sterile, 30 kits/cs)

... includes a catheter adapter, an injection cap and a clamp.

QC Cap QC-200 (Prefilled with Povidone Iodine, individually packed, 30 pkgs/cs)

... used to cap off the catheter adapter. It provides a bacteriological safeguard for CAPD patients in between exchanges.

Felt Cuff FC-1000 (Sterile, 1 cm long, 10 cuffs/package)

... enables the physician to custom fit peritoneal catheters. It is designed to fit over silicone tubes with nominal outer diameter of 4.9 mm.

Medical Grade Silicone Adhesive MGS-520 (Type A, 8 gm, sterile, 1/package)

... is used to bond felt cuffs to a catheter. It is also used to repair or extend peritoneal catheters.

Peritoneal Catheter Repair Kit RC-1000 (Sterile, 1/package)

... is used to replace or repair the external portion of a damaged catheter. It contains -
1-6" long silicone catheter tubing with sleeve and adapter
1-catheter clamp #PCA-1030
1-injection cap #ICL-100

Titanium Adapter TA-200 (Sterile, 1/package)

... is a two part titanium connector with a locking nut. It is used to terminate the external portion of peritoneal catheters. It allows peritoneal catheters to connect securely to CAPD administration sets.

Adult Titanium Catheter Extender TE-1800 (Sterile with 6" long silicone tubing, 1/package)

... includes 6" long silicone catheter tubing with a Titanium Extender. It is used to repair or extend the external portion of adult peritoneal catheters.

Pediatric Titanium Catheter Extender TE-1800P (Sterile with 6" long silicone tubing, 1/package)

... includes 6" long silicone catheter tubing with a Titanium Extender. It is used to repair or extend the external portion of peritoneal catheters with pediatric sized tubing.

Faller Stylet (Trocar) FS-401 (Sterile, autoclavable, 1/package)

... used to create subcutaneous tunnel. It contains a cutting edge on one end and a barbed section on the other end to connect to the catheter. Designed to fit Tenckhoff catheters.

Trocar TR-400 (Non-sterile, autoclavable, 1/package)

... used to insert PD catheters. It is a tapered bivalve instrument with a handle on each half and an obturator with a trocar point. Designed to fit standard Tenckhoff or Reinforced Tenckhoff PD catheters.

Stylet SST-xxx (Non-sterile, autoclavable, 1/package)

... is used to straighten the peritoneal catheters. It is available in different lengths for use with catheters of different overall lengths.
Standard lengths are -
SST-400 (400 mm), SST-445 (445 mm), SST-470 (470 mm), SST-500 (500 mm), SST-650 (650 mm)

Peritoneal Tunneling Stylet PS-423 (Sterile, single use only, 1/package)

... is made of pliable stainless steel which can be bent to create a customized curved subcutaneous tunnel.

16 Fr Pull-Apart Introducer Set PS-16FR (Sterile, 1/package)

... is used to insert catheters percutaneously using modified Seldinger technique. It contains one 16 Fr Pull-Apart Sheath/Dilator, one 18ga Introducer Needle, one 10cc Syringe and one J/Straight 0.038" Guidewire.